



**ELITE/SAVVY 8000 Series ® Air Mattress System
User Manual**



Air Kinetic Technologies Corp.

TEL : +886-3-4194567

FAX : 886-3-4693789

Address : No. 3, Gongye 5thRd., Pingzhen 324, Taoyuan, Taiwan

<http://www.aktc.com.tw>

IMPORTANT SAFEGUARDS

READ ALL INSTRUCTIONS BEFORE OPERATING THIS DEVICE



NOTE, CAUTION AND WARNING STATEMENTS:

NOTE –Indicate some tips

CAUTION – Indicate correct operating or maintenance procedures in order to prevent damage to or destruction of the equipment or other property

WARNING – Call attention to a potential danger that requires correct procedures or practices in order to prevent personal injury.



WARNING – To reduce the risk of electrocution

1. Always unplug this product immediately while it's not in use.
2. Do not disassemble the pump to avoid electrocution
3. Do not place or store product where it can fall or be pulled into a tub or sink.
4. Do not place in or drop into water or other liquid. Do not use while bathing.
5. Do not reach for a product that has fallen into water. Unplug immediately.



WARNING – To reduce the risk of burns, electrocution, fire or injury to persons

1. The operation of the system has to have the mattress connected to the PUMP, please do not power-off or unplug the PUMP in operation.
2. This product should never be left unattended when plugged in.
3. Close supervision is necessary when this product is used by, on, or near children or invalids.
4. Use this product only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
5. Never operate this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Return the product to a service center or to the distributor for examination and repair.
6. Keep the cord away from heated surfaces.
7. Never block the air openings of this product or place it on a soft surface, such as a bed or couch, where the openings may be blocked. Keep the air opening free of lint, hair, and other similar particles.
8. Never drop or insert an object into any opening or hose.
9. Connect this product to a properly grounded outlet only. See Grounding Instruction.
10. Put the power cord or hose tube at the patient foot area to avoid wound on the patient's head.
11. To avoid electromagnetic interference, the patient environment should not have










strong electro-magnetic or RF generated equipment near by

12.The PUMP will have minor heat generated in operation, please do not direct contact the surface continuously for more than 1 minute.

13.The product with ground pin (3 pin type) cannot be used in the home.

14.The EMC specification is compliant with the regulation requirement (please reference to the EMC information at the last page). For power cord with ground pin (3 pin type), the connection with properly grounded power outlet would get a better EMC suppressing effect. The system will work correctly also for the power cord connection with the power outlet without grounding.

15.When loss or failure of the supply mainstemporarily(In 20 minutes). It causespump stop, fail of power indication and alarm. But these are normal. The product can return to work state after supply mains is stable.

SYMBOLS	DESCRIPTION
I	POWER ON
O	POWER OFF
	ATTENTION
	DOUBLE ISOLATION
	“BF” SYMBOL, INDICATE THIS PRODUCT IS ACCORDING TO THE DEGREE OF PROTECTION AGAINST ELECTRIC SHOCK FOR TYPE BF EQUIPMENT
	CAUTION, READ THE INSTRUTION BEFORE USE
	AWAY FROM THE FLAME
IP21	WATER AND DUST PROTECTION CLASSIFICAITON
	FUSE SPECIFICATION
	DISPOSAL OF ELECTRICAL & ELECTRONIC EQUIPMENT(WEEE): THIS PRODUCT SHOULD BE HANDED OVER TO AN APPLICABLE COLLECTION POINT FOR THE RECYCLING OF ELECTRICAL AND ELECTRONIC EQUIPMENT.
	UL CERTIFICATION LOGO (COMPLIACE WITH IEC60601-1) With respect to electrical shock, fire and mechanical hazards only in accordance with STANDARD.
	CB CERTIFICATION LOGO



Contents

1 INTRODUCTION 5

2 INTENDED USE..... 5

3 PRODUCT DESCRIPTION 5

4 PRODUCT INSTALLATION GUIDE..... 5

5 PANEL DISPLAY AND THE OPERATION GUIDE 8

6 CLEANING 13

7 STORAGE..... 14

8 MAINTENANCE 14

9 The Disposal of Air Mattress 14

10 TROUBLESHOOTING 15

11 TECHNICAL DATA 16

1 INTRODUCTION

This manual should be used for initial set up of the AKTC **ELITE/SAVVY 8000 Series® Air Mattress System** and for daily maintenance. **Please keep the manual in handy area for reference.**

2 INTENDED USE

This product is intended to help and reduce the incidence of pressure ulcers while optimizing patient comfort. It also provide following purposes:

- Individual home care setting and long-term care of whom suffering from pressure ulcer.
- Pain management as prescribed by physician.

⚠ NOTE: Equipment not suitable for use in the presence of a flammable anesthetic mixture with air or with oxygen or nitrous oxide.

3 PRODUCT DESCRIPTION


The ELITE/SAVVY 8000 Series are alternating mattress replacement system used in the prevention and treatment of pressure ulcers. By using the established principles of alternating therapy, the ELITE/SAVVY 8000 Series offer patients a comfortable and relaxing support surface which can both prevent skin breakdown and enhance healing.

The CONTROL UNIT of the ELITE/SAVVY 8000 Series is a compact pump featuring an audible and visual low pressure, power failure and machine malfunction alarms, and a digital pressure adjustment function. The 19 cells mattress unit provides a unique design which keeps the lower layer of air cells constantly inflated while alternating and deflating the upper layer. The head section of cells remains static. The mattress has a heavy-duty nylon base sheet with a vapor permeable PU coated stretch cover.







In the event of cardiac arrest, rapid deflation is achieved by using the highly visible CPR facility.

4 PRODUCT INSTALLATION GUIDE

1. Unpacking the box to inspect for any damage, which may have occurred during shipment. If there are any damages, please contact your dealer immediately.
2. Place the mattress on top of the bed frame. Please note for the foot end.

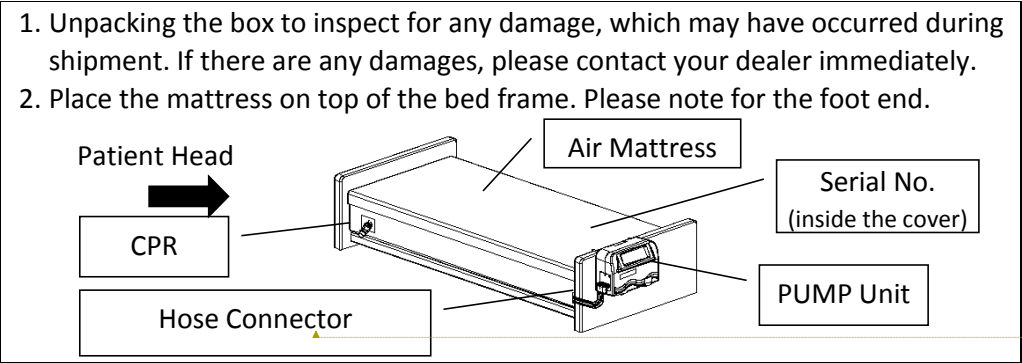


Formatted: Font: 14 pt, Complex Script
Font: +Body (Calibri), 14 pt
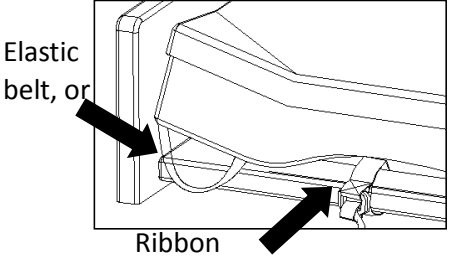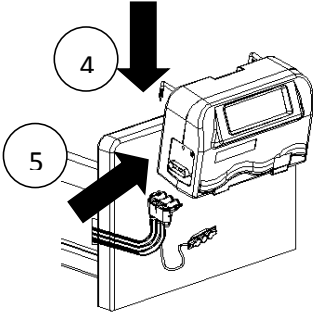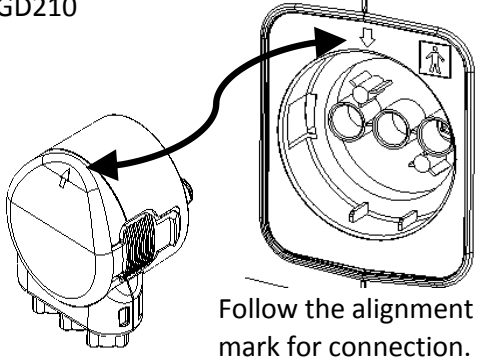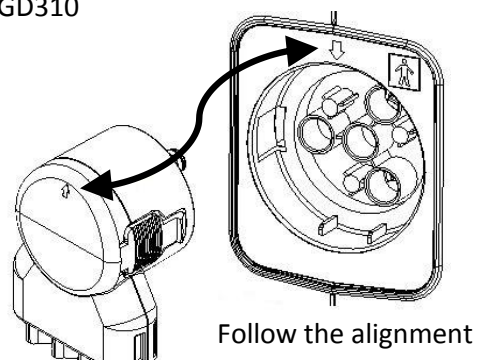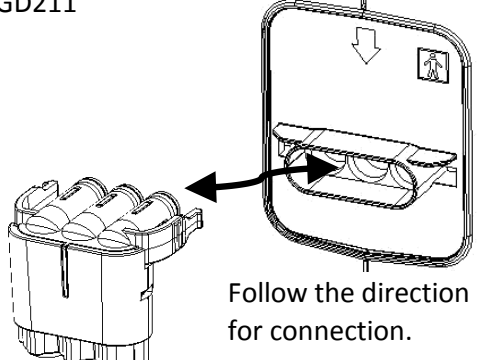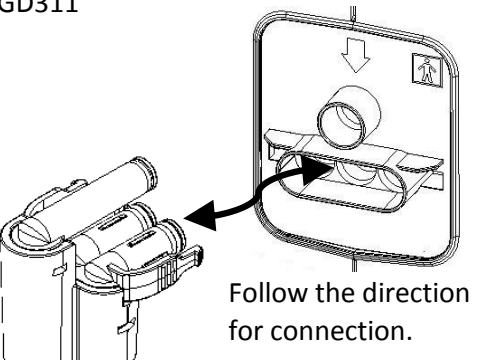
<p>3. Fix the mattress onto the bed frame</p>	 <p>Elastic belt, or Ribbon</p>
<p>4. Hang the pump onto bed rail (foot-end), the hangers will hold the bed rail tight automatically</p> <p>5. Unplug the cover of the hose connector and connect the hose connector to the pump unit.</p> <p>Note: There' s 4 combinations of hose connectors vs. the pumps' side panel as following (GD210, GD211, GD310, GD311), make sure you hear the click sound at connection.</p>	
<p>GD210</p>  <p>Follow the alignment mark for connection.</p>	<p>GD310</p>  <p>Follow the alignment mark for connection.</p>
<p>GD211</p>  <p>Follow the direction for connection.</p>	<p>GD311</p>  <p>Follow the direction for connection.</p>


6. Plug the power cord into electrical outlet

- ⚠ NOTE: Make sure the pump unit is suitable for the local power voltage
- ⚠ CAUTION: The pump can only be applied to the mattress recommended by the manufacturer. Do not use it for any other purpose (applied part: air mattress)

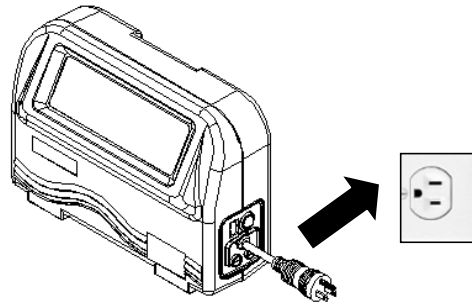


7. Make sure the CPR is at CLOSE position before turning on the power. Switch the CPR to OPEN position to release the air at emergency or for packaging

CPR Position Indicator

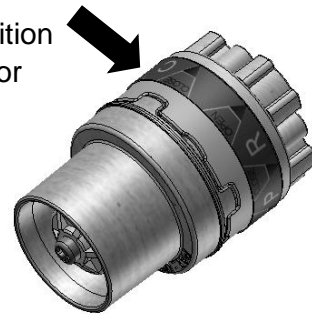


8. Put on the hose connector cover at transportation, the mattress will retain pressure for up to 24hrs


GD210
After

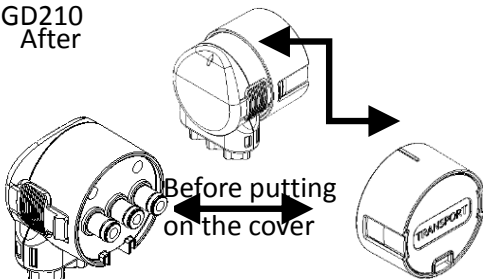

GD310

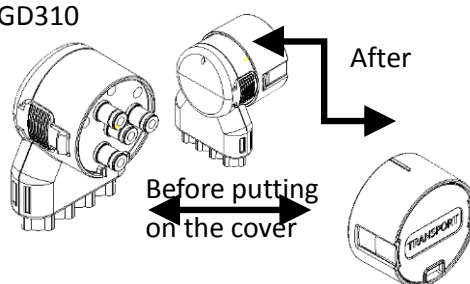


GD211

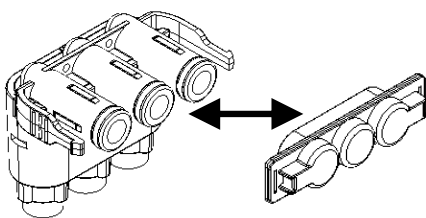


No direction for the cover

GD311

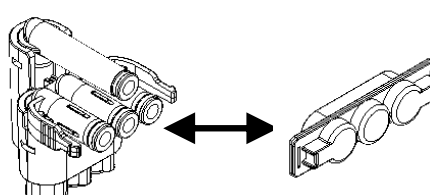

No direction for the cover

5 PANEL DISPLAY AND THE OPERATION GUIDE

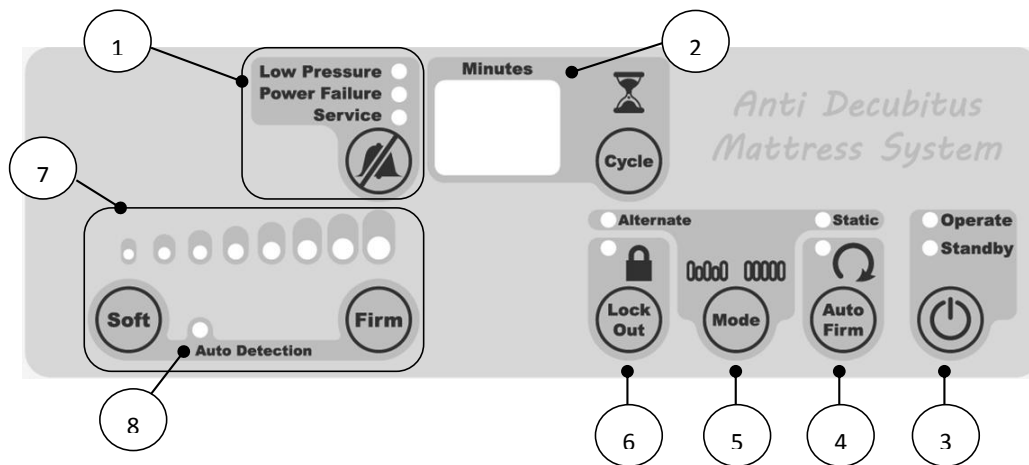

5.1 PANEL DISPLAY

- ① Alarm Mute and Alarm Indicator
 - Low Pressure Alarm Indicator
 - Power Failure Alarm Indicator
 - Service (Malfunction) Alarm Indicator
- ② Alternate Cycle Time or Warning code Display
- ③ Operating or Standby
- ④ Auto-Firm
- ⑤ Function Mode Selection (Alternate & Static)
- ⑥ Panel Lock-out
- ⑦ Comfort Control
- ⑧ Auto Detection



5.1.1 ALARM MUTE

Press alarm mute button to temporary suspend the Low-Pressure/Power Failure/Static Overtime /Service alarms. Should the situation not resolved and the fault conditions continue, the alarm shall resume notifying the patient/caregiver.



5.1.2 Alternate Cycle Time Display

Alternating Cycle Time can be selected from 10~30mins at 5mins interval by pressing the CYCLE button



5.1.3 Operate or Standby

Press this button to start operating or go into standby.

NOTE: The power switch on the side of pump must be turned on. At Power on the unit will resume the state before last power-off.



5.1.4 Inflate/Auto-Firm

The PUMP will go into the Inflate mode (LED lights flashing) every time the OPERATE mode is triggered. This insures the mattress to be able to reach its maximum operating pressure. Once the max pressure level is reached, the pump will automatically switch into the previous selected mode and comfort level. User can also use this function as full mattress inflation during patient sit-up or ingress/egress for better support.



5.1.5 Function Mode Switch

- **ALTERNATE** - for the mattress to operate at alternating mode the air cell of the mattress will be proportionally deflated to reduce the surface pressure. The alternating cycle will continue at the selected cycle time until another mode is selected.
- **STATIC** - This mode allows the mattress to maintain at the selected pressure. After 20 minutes, the **STATIC OVER TIME** alarm will be triggered for 10mins at every 15 seconds interval. Without further action, the PUMP will go into **ALTERNATE** mode automatically.



5.1.6 Panel Lock-Out


Should the panel remain untouched for 30 seconds or press the Lock-out button, the lock-out feature will lock the screen to prevent accident from changing the setting without notice. To unlock, press the Lock-out button for 3 seconds.



5.1.7 Comfort Level

Comfort level controls the air pressure output. When pressing the **FIRM** button,

the output pressure will increase and higher pressure output will support the heavier weight user, for decreasing air pressure, vice versa. Check to see if the suitable pressure is selected by sliding one hand between the air cells and the patient to feel patient's buttocks. Users should be able to feel the minimum contact. Always leave at least 1 inch space between user's buttock areas and air cells under to prevent bottoming out.

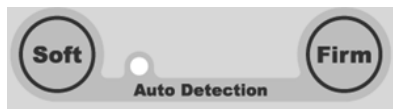

5.1.8 Auto Detection


When pressing the SOFT and FIRM button together, the pump will automatically detect the weight of the patient and set the appropriate pressure output for patient comfort.


5.2 OPERATION GUIDE

5.2.1 GENERAL OPERATION :



NOTE: The power switch is located on the side of pump

- Press  to turn on the unit, all indicators on the control panel will light up accompanied with a beep for 2 seconds (You can also check the indicator for failure if any), and the indicator of STANDBY on the control panel will light up (In case the pump was turning off at OPERATE, it will go to OPERATE directly).

Ps: To test if the battery is working properly, press  to turn off the power. Power failure alarm should be triggered. If not, please call customer service.












- Push on the OPERATE button , the system will start inflation and the "AUTO-FIRM" indicator will be flashing.
- The mattress should be fully inflated within 60 minutes, and automatically enter the last operating mode, otherwise the low pressure alarm will be triggered.
- According to the weight of the patient, adjust the pressure setting to the most suitable level without bottoming out. User can determine an appropriate pressure by adjusting the Comfort Level. Please consult with your physician for a proper setting.

Table 1 Weight and Comfort Level Reference Table


GDseries pump+8" mattress (iLAL/2-1Alternate)										
Comfort Control (Auto-Detection)	Pump output Pressure(mmHg)	Patient Weight (KG)								
		20	40	60	80	100	120	140	160	180
	25	< 40								
	30		20~60							
	35			40~80						
	40				60~100					
	45					80~120				
	50						100~140			
	55							120~160		
	60								140~180	

5.2.2 CPR

When CPR needs to be performed, quickly rotate the CPR valve to "OPEN" position, at the same time, disconnect the hose connector from the PUMP to speed up the air release.

5.2.3 AUDIBLE AND VISIBLE ALARM






- (1) Power Failure – When electrical shortage occurred or power cord is unplug without turning off the pump, the "POWER FAILURE" indicator will light up along with buzzer. Check to ensure power cord is connected properly

 **NOTE:** When the PUMP is not been used for more than 3 months, it might need 6hours operating time or more for the Alarm to function properly.

- (2) Low Pressure –When an abnormal low pressure occurred in body section for 4mins, the "Low Pressure" indicator will flash and beep every 4 seconds. The Low Pressure alarm will continue until alarm mute button being pressed. Should the situation not be resolved and the fault conditions continue, the alarm will resume.
- (3) Static Overtime –When Static mode lasts for more than 20 minutes, "Static" indicator will flash and beep every 15 seconds. Press the MODE key to change to alternating mode can disable the alarm, or press the MUTE button to pause the alarm for 20 minutes. If not dealt, the system will automatically enter into alternating mode after 10 minutes.
- (4) Power Failure Overtime – When the power is lost for more than 20mins, the "Power Failure" indicator will flash along with a beep sound at every 15 seconds, Press the MUTE button to cease the alarm, and check the patient for bed sore examination.

- (5) Service(Malfunction)–When fault conditions occur, the "SERVICE" indicator will light up along with buzzer. Reference to Table 2 for Warning code and call the agent or distributor for service.

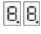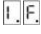Table 2 WarningCode Reference Table

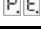PRIORITY HIHG ↓ LOW	WARNINGCODE	INDICATOR LED	AUDIBLE OUTPUT MODE	CONDITION OF OUTPUT	WARNING DESCRIPTION	REMARKS
0	N/A	N/A	ONCE	Not in System Shutdown	Key Tone	Key Tone from Functional Button
1		Power Failure	ONCE	POWER-OFF	System Shutdown	
2		ALL LED	ONCE	OPERATE OR STANDBY	Power-On	All Indicators Light On
3	N/A	N/A	ONCE	OPERATE OR STANDBY	State/Mode Switching	
4		AutoFirm	ONCE	OPERATE	Mattress Inflation Completion	Inflation Ended
5		AutoFirm	ONCE	OPERATE	Auto-Firm Completion	Auto-Firm Ended
6		Static	ONCE	OPERATE	Static Completion	Static Ended
7	N/A	Power Failure	REPEAT (cycle4sec.)	POWER-OFF	Power Failure Alarm	No Display
8		Low Pressure	REPEAT (cycle4sec.)	OPERATE OR STANDBY	Power-On Inflation FailureAlarm	
9		Low Pressure	REPEAT (cycle4sec.)	OPERATE OR STANDBY	Auto-Firm Failure Alarm	
10		Low Pressure	REPEAT (cycle4sec.)	OPERATE OR STANDBY	Low Pressure Overtime Alarm	
11		Service	REPEAT (cycle4.5sec.)	OPERATE OR STANDBY	Constant Pressure Control Failure Alarm	
12		Service	REPEAT (cycle4.5sec.)	OPERATE OR STANDBY	High Pressure Overtime Alarm	
13		Service	REPEAT (cycle4.5sec.)	OPERATE OR STANDBY	Low Ambient Temperature Alarm	Environment Temperature Over Specification Limit
14		Service	REPEAT (cycle4.5sec.)	OPERATE OR STANDBY	High Ambient Temperature Alarm	Environment Temperature Over Specification Limit
15		Service	REPEAT (cycle4.5sec.)	OPERATE OR STANDBY	Air Valve 1 Positioning Failure Alarm	Air Valve 1 failure
16		Service	REPEAT (cycle4.5sec.)	OPERATE OR STANDBY	Air Valve 2 Positioning Failure Alarm	Air Valve 2 failure
17		Power Failure	REPEAT (cycle15sec.)	OPERATE OR STANDBY	Power Failure Overtime Alarm	
18		Service	REPEAT (cycle15sec.)	OPERATE OR STANDBY	Battery Low Alarm	Battery would need to be replaced
19		Static	REPEAT (cycle15sec.)	OPERATE OR STANDBY	Static Overtime Alarm	
20		NONE	NONE	FACTORY CALIBRATION MODE	Calibration Not Completed	
21		NONE	NONE	FACTORY CALIBRATION MODE	Calibration Completed	


5.2.4 ALARM MUTE

When alarms were triggered, both the LED light and buzzer will sound off to warn the patient/caregiver. By pressing the button, it will temporary mute the buzzer so the caregiver may check for possible causes. Should the situation not resolved

and fault conditions continue, the alarm will resume. When in Power Failure situation, pressing alarm mute will cease the buzzer and turn off the "Power Failure" indicator.

6 CLEANING

By wiping the PUMP UNIT with a damp cloth pre-soaked with a mild detergent, and keep it away from dust. If other detergent is used, choose one that will have no chemical effects on the surface of the plastics case of the pump unit.

 **CAUTION:** Do not immerse or soak pump unit.

By using a single use wipe, clean the MATTRESS COVER with a solution of neutral detergent and hand hot water. Rinse thoroughly with clean water and a damp single use wipe.

Disinfecting the cover

If the cover is heavily soiled or has been exposed to bodily fluids such as blood, it will require a more thorough cleaning procedure.

Wipe the cover using a single use wipe and a 0.1% Chlorine Solution (1,000ppm) and cold water. If required a 1% Chlorine Solution (10,000ppm) and cold water can be used. Rinse thoroughly with clean water and a damp single use wipe. Make sure the cover is completely dried before refitting to the mattress.

Frequent or prolonged exposure to higher concentration disinfectant solutions may prematurely age the fabric cover of mattresses. Surfaces must be protected during use and rinsed and thoroughly dried after application of a disinfectant.

Laundering

Before laundering mattress covers should be completely removed. Where required mattress covers can be laundered as follows:

Pre wash 60°C + 15 minutes

Main wash 60°C + 15 minutes


This should be followed by a cold rinse and extraction.

Drying

Mattress covers should be hung from a line or bar and drip dried in a clean indoor environment. Covers must be completely dried before refitting to the mattress.

Mattress covers can be tumble dried on a low heat setting for 90 minutes. Drying

temperature must not exceed 40°C . Exceeding the temperature can cause significant damage to the mattress cover.

 **CAUTION:** Do not use phenolic-based product for cleaning.

 **CAUTION:** After cleaning, dry the mattress without direct exposure of sunlight.

7 STORAGE

- To quickly vacuum air out from mattress for storage, rotate the CPR valve to OPEN position and disconnect the hose connector to release the air.
- Lay the mattress out flat and upside down.
- Roll from the head end towards the foot end
- Foot-end strap can then be stretched around the rolled mattress to prevent unrolling
- The power cord could be wrapped around the pump bumper or disconnected for storage.

8 MAINTENANCE

8.1 8.1 General

- Check main power cord and plug if there are abrasions or excessive wears.
- Check mattress cover for signs of wear or damage. Ensure mattress cover and tubes are stubbed together correctly.
- Check the air hoses for any kink or break. For replacement, please contact your local dealers.

8.2 8.2 FUSE REPLACEMENT

- Disconnect the plug from mains power when a blown fuse is suspected.
- Remove the cover of the fuse holder by means of a small screwdriver.
- Insert a new fuse of the correct rating in, and replace the cover of the fuse holder back. The fuse rating should comply with the requested specification.

8.3 8.3 AIR FILTER REPLACEMENT

- Replace the air filter located at the back of the pump.
- The filter is reusable and can be washed gently with a mild detergent and water. Dry the filter before use.
- Check and replace air filter regularly if environment is dirty.

9 The Disposal of Air Mattress

When the air mattress is broken or no longer be useable, the mattress and the pump may be discard for recycle.


10 TROUBLESHOOTING

PROBLEM	SOLUTION
The mattress is not able to connect with the PUMP	<ul style="list-style-type: none"> ● Check if the mattress model (model no. located inside the cover at the foot end)xxAAAxix matched with the PUMP model xxBBB-xxx. The AAA should be the same as BBB. If not, please contact with the agent or distributor ● Check if the connector cover is removed and make sure the connector is not broken
The pump is showing no indications it is working	<ul style="list-style-type: none"> ● Check if the plug is connected to mains ● Check if the main power switch is at ON position ● Check if there is any blown fuse
Power Failure Alarm Failure	<ul style="list-style-type: none"> ● The pump is in operation but the power failure alarm is not working at power down, please call customer service
The low pressure light is constantly flashing and the alarm is sounded	<ul style="list-style-type: none"> ● Check if the CPR is at CLOSE position ● Check if the power was suddenly shut down ● Check if the connection between air tube to pump unit is tightly secured ● Check if all coupling connections along mattress are secured
Power Failure Alarm Failed	<ul style="list-style-type: none"> ● If the PUMP is in operation but failed to trigger the Power Failure Alarm at Power Off, please contact the dealer or agent for further investigation
The pump is on but the mattress is not alternated	<ul style="list-style-type: none"> ● Make sure the mattress inflationis completed ● Check the pump control panel the indicator of “ALTERNATE” is lighted on, if not, switch it to “ALTERNATE” ● Check if “Service” alarm indicator is on with buzzer, if yes, contact the dealer or agent for further investigation
The pump is operating noisily	<ul style="list-style-type: none"> ● Make sure the pump is resting against a solid surface ● If the noise getting louder, contact the dealer or agent for further investigation
Patient is bottoming out (without alarm triggered)	Pressure setting might be inadequate for the patient, adjust comfort level to FIRM and wait for a few minutes for better comfort

If the above information does not solve the problem, please contact your localdealer or agent for further support.

11 TECHNICAL DATA

11.1 Product Specification

PUMP UNIT		AIR MATTRESS	
MODEL	GD Series Pump	MODEL	8" Mattress Series
DIMENSION(cm)	33 (W) x 22 (D) x12 (H)	DIMENSION(cm)	89 (W) x 200 (L) x 21 (H)
WEIGHT(kg)	3.5kg	WEIGHT(kg)	10Kg
CYCLE TIME	10/15/20/25/30minutes	CELL MATERIAL	Nylon TPU or TPU film
STATIC TIME	30 minutes		
AUTO FIRM TIME	20 minutes		
PUMP OUTPUT FLOW RANGE (Liter)	> 8L (@120V or 230V) Note: The flow rate may be varied because of the fluctuation of input voltage		
PUMP OUTPUT PRESSURE RANGE (mmHg)	25 to 60 (± 2)	NO. OF AIR CELL	19 CELLS
POWER	AC120V 60Hz or AC230 V 50Hz	COVER MATERIAL	2Way Stretched Polyester with PU coated
CURRENT	0.25 A _{MAX} (@132V~) or 0.12A _{MAX} (@253V~)	BOTTOM MATERIAL	Nylon-TPU
FUSE RATING	T1AL 250VAC	MAX WEIGHT	180kg
FREQUENCY	60Hz (120V) or 50 Hz (230V)	MAX PRESSURE	103.5 mmHg
CLASSIFICATION	Class II		
	Type BF 		
WARRANTY	2 years	WARRANTY	1 year
SHELFLIFE	2 years	SHELFLIFE	1 year
ENVIRONMENTAL CONDITIONS			
OPERATION ENVIRONMENT	5°C~40°C 15%RH ~ 93%RH(no condensation)		
STORAGE ENVIRONMENT	-25°C~70°C ≤93%RH(no condensation)		
ENVIRONMENT PRESSURE	70 kPa-101.3kPa		
ENVIRONMENTHORIZONTAL LEVEL	≤ 3000m		
WATER AND DUST PROTECTION CLASSIFICAITON	IP21		

11.2 EMC INFORMATION

Guidance and manufacturer’s declaration-electromagnetic emissions		
The GD311-301 is intended for use in the electromagnetic environment specified below. The customer or the user of the GD311-301 should assure that it is used in such an environment.		
Emission test	Compliance	Electromagnetic environment-guidance
RF emissions CISPR 11	Group 1	The GD311-301 uses RF energy only for its internal function. Therefore, its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment. The GD311-301 is suitable for use in all establishments, including domestic establishments and those directly connected to the public low-voltage power supply network that supplies buildings used for domestic purposes.
RF emissions CISPR 11	Class B	
Harmonic emissions IEC 61000-3-2	Class A	
Voltage fluctuations /flicker emissions IEC 61000-3-3	Compliance	


Guidance and manufacturer's declaration-electromagnetic immunity

The [GD311-301](#) is intended for use in the electromagnetic environment specified below.

The customer or the user of the [GD311-301](#) should assure that it is used in such an environment.

Immunity test	IEC 60601 test level	Compliance level	Electromagnetic environment-guidance
Electrostatic discharge(ESD) IEC 61000-4-2	± 6 kV contact ± 8 kV air	± 6 kV contact ± 8 kV air	Floors should be wood, concrete or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30%
Electrical fast transient/burst IEC 61000-4-4	± 2kV for power supply lines ± 1kV for input/output lines	± 2kV for power supply lines Not applicable	Mains power quality should be that of a typical commercial or hospital environment.
Surge IEC 61000-4-5	± 1kV line(s) to line(s) ± 2kV line(s) to earth	± 1kV differential mode Not applicable	Mains power quality should be that of a typical commercial or hospital environment.
Voltage Dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11	<5% UT(>95% dip in UT) for 0,5 cycle 40% UT(60% dip in UT) for 5 cycles 70% UT(30% dip in UT) for 25 cycles <5% UT(>95% dip in UT) for 5 s	<5% UT(>95% dip in UT) for 0,5 cycle 40% UT(60% dip in UT) for 5 cycles 70% UT(30% dip in UT) for 25 cycles <5% UT(>95% dip in UT) for 5 s	Mains power quality should be that of a typical commercial or hospital environment. If the user of the GD311-301 requires continued operation during power mains interruptions, it is recommended that the GD311-301 be powered from an uninterruptible power supply or a battery.
Power frequency(50/60 Hz) magnetic field IEC 61000-4-8	3 A/m	3 A/m	The GD311-301 power frequency magnetic fields should be at levels characteristic of a typical location in a typical commercial or hospital environment.

NOTE UT is the a.c. mains voltage prior to application of the test level.

Guidance and manufacturer's declaration-electromagnetic immunity			
The GD311-301 is intended for use in the electromagnetic environment specified below.			
The customer or the user of the GD311-301 should assure that is used in such and environment.			
Immunity test	IEC 60601 test level	Compliance level	Electromagnetic environment-guidance
Conducted RF IEC 61000-4-6	3 Vrms 150 KHz to 80 MHz	3 Vrms	<p>Portable and mobile RF communications equipment should be used no closer to any part of the GD311-301 including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter.</p> <p>Recommended separation distance:</p> $d = 1,2 \sqrt{P}$ $d = 1,2 \sqrt{P} \quad 80\text{MHz to } 800 \text{ MHz}$ $d = 2,3 \sqrt{P} \quad 800\text{MHz to } 2,5 \text{ GHz}$ <p>Where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and d is the recommended separation distance in metres (m).</p> <p>Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey,^a should be less than the compliance level in each frequency range.^b</p> <p>Interference may occur in the vicinity of equipment marked with the following symbol:</p> 
NOTE1 At 80 MHz and 800 MHz, the higher frequency range applies. NOTE2 These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.			
<p>a Field strengths from fixed transmitters, such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location in which the GD311-301 is used exceeds the applicable RF compliance level above, the GD311-301 should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as re-orienting or relocating the GD311-301.</p> <p>b Over the frequency range 150 kHz to 80 MHz, field strengths should be less than 3 V/m.</p>			

Recommended separation distance between portable and mobile RF communications equipment and the GD311-301			
The GD311-301 is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. The customer or the user of the GD311-301 can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and the GD311-301 as recommended below, according to the maximum output power of the communications equipment.			
Rated maximum output power of transmitter W	Separation distance according to frequency of transmitter m		
	150 kHz to 80 MHz $d = 1,2\sqrt{P}$	80 MHz to 800 MHz $d = 1,2\sqrt{P}$	800 MHz to 2,5 GHz $d = 2,3\sqrt{P}$
0,01	0,12	0,12	0,23
0,1	0,38	0,38	0,73
1	1,2	1,2	2,3
10	3,8	3,8	7,3
100	12	12	23
For transmitters rated at a maximum output power not listed above, the recommended separation distance d in metres (m) can be estimated using the equation applicable to the frequency of the transmitter, where p is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer.			
NOTE1 At 80 MHz and 800 MHz, the separation distance for the higher frequency range applies.			
NOTE2 These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.			


Guidance and manufacturer's declaration-electromagnetic emissions		
The GD311-401 is intended for use in the electromagnetic environment specified below.		
The customer or the user of the GD311-401 should assure that it is used in such an environment.		
Emission test	Compliance	Electromagnetic environment-guidance
RF emissions CISPR 11	Group 1	The GD311-401 uses RF energy only for its internal function. Therefore, its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment.
RF emissions CISPR 11	Class B	
Harmonic emissions IEC 61000-3-2	Class A	
Voltage fluctuations /flicker emissions IEC 61000-3-3	Compliance	

Guidance and manufacturer's declaration-electromagnetic immunity

The [GD311-401](#) is intended for use in the electromagnetic environment specified below.

The customer or the user of the [GD311-401](#) should assure that it is used in such an environment.

Immunity test	IEC 60601 test level	Compliance level	Electromagnetic environment-guidance
Electrostatic discharge(ESD) IEC 61000-4-2	± 6 kV contact ± 8 kV air	± 6 kV contact ± 8 kV air	Floors should be wood, concrete or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30%
Electrical fast transient/burst IEC 61000-4-4	± 2 kV for power supply lines ± 1 kV for input/output lines	± 2 kV for power supply lines Not applicable	Mains power quality should be that of a typical commercial or hospital environment.
Surge IEC 61000-4-5	± 1 kV line(s) to line(s) ± 2 kV line(s) to earth	± 1 kV differential mode Not applicable	Mains power quality should be that of a typical commercial or hospital environment.
Voltage Dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11	<5% UT(>95% dip in UT) for 0,5 cycle 40% UT(60% dip in UT) for 5 cycles 70% UT(30% dip in UT) for 25 cycles <5% UT(>95% dip in UT) for 5 s	<5% UT(>95% dip in UT) for 0,5 cycle 40% UT(60% dip in UT) for 5 cycles 70% UT(30% dip in UT) for 25 cycles <5% UT(>95% dip in UT) for 5 s	Mains power quality should be that of a typical commercial or hospital environment. If the user of the GD311-401 requires continued operation during power mains interruptions, it is recommended that the GD311-401 be powered from an uninterruptible power supply or a battery.
Power frequency(50/60 Hz) magnetic field IEC 61000-4-8	3 A/m	3 A/m	The GD311-401 power frequency magnetic fields should be at levels characteristic of a typical location in a typical commercial or hospital environment.
NOTE UT is the a.c. mains voltage prior to application of the test level.			

Guidance and manufacturer's declaration-electromagnetic immunity			
The GD311-401 is intended for use in the electromagnetic environment specified below.			
The customer or the user of the GD311-401 should assure that is used in such and environment.			
Immunity test	IEC 60601 test level	Compliance level	Electromagnetic environment-guidance
Conducted RF IEC 61000-4-6	3 Vrms 150 KHz to 80 MHz	3 Vrms	<p>Portable and mobile RF communications equipment should be used no closer to any part of the GD311-401 including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter.</p> <p>Recommended separation distance:</p> $d = 1,2 \sqrt{P}$ $d = 1,2 \sqrt{P} \quad 80\text{MHz to } 800 \text{ MHz}$ $d = 2,3 \sqrt{P} \quad 800\text{MHz to } 2,5 \text{ GHz}$ <p>Where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and d is the recommended separation distance in metres (m).</p> <p>Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey,^a should be less than the compliance level in each frequency range.^b</p> <p>Interference may occur in the vicinity of equipment marked with the following symbol:</p> 
NOTE1 At 80 MHz and 800 MHz, the higher frequency range applies.			
NOTE2 These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.			
a Field strengths from fixed transmitters, such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location in which the GD311-401 is used exceeds the applicable RF compliance level above, the GD311-401 should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as re-orienting or relocating the GD311-401 .			
b Over the frequency range 150 kHz to 80 MHz, field strengths should be less than 3 V/m.			

Recommended separation distance between portable and mobile RF communications equipment and the GD311-401			
The GD311-401 is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. The customer or the user of the GD311-401 can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and the GD311-401 as recommended below, according to the maximum output power of the communications equipment.			
Rated maximum output power of transmitter W	Separation distance according to frequency of transmitter m		
	150 kHz to 80 MHz $d = 1,2\sqrt{P}$	80 MHz to 800 MHz $d = 1,2\sqrt{P}$	800 MHz to 2,5 GHz $d = 2,3\sqrt{P}$
0,01	0,12	0,12	0,23
0,1	0,38	0,38	0,73
1	1,2	1,2	2,3
10	3,8	3,8	7,3
100	12	12	23
For transmitters rated at a maximum output power not listed above, the recommended separation distance d in metres (m) can be estimated using the equation applicable to the frequency of the transmitter, where p is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer.			
NOTE1 At 80 MHz and 800 MHz, the separation distance for the higher frequency range applies.			
NOTE2 These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.			